Aidite®

New Intraoral Scanner

cameo Elegant



Meet cameo Elegant 3

Empowering a Highly Efficient Digital Chairside Solution

Cameo Elegant 3

Say hello to cameo brand new intraoral scanner, Elegant 3! Exceptional, smart and elegant, Elegant 3 improves the overall clinical experience. It allows users to easily acquire vivid digital impressions in a faster and more accurate way.

Moreover, Elegant 3 comes with a wide range of different intelligent functions, guaranteeing more comfortable chairside experience and efficient clinic-lab collaboration.

Ergonomically Designed for Optimal Hand-held Operation

Dynamic LED Indicator

Automatic Anti-fogging


Autoclavable(100 times)

One-button Control Montion Sensing Control

Light Weight: 240g





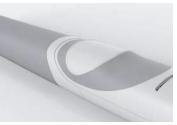





Capture Digital Impressions in Real Time

Fast, Accurate and Exceptional Scan Performance for Chairside Experience

The Elegant 3 scanner optimizes the hardware structure and software algorithm to make scanning faster.



Faster Scan

The scan speed is 30% faster than the last generation due to Elegant 3 drastically upgraded scanning technology.



More Accurate Scan

Thanks to optimized algorithms and imaging mechanism, Elegant 3 carries out more accurate and precise results by 30%, meeting the clinical needs of different applications.



Longer Scan Tip

Slimmer and 15% longer autoclavable scanner tips offer patients with more comfortable clinical experience.



Greater Scan Depth

Elegant 3 is suitable for various clinical scenarios including acquiring scan bodies and periodontal scans.



User-Friendly


Smart and Intuitive - Smooth Scanning Experience for Dentists

With powerful and intelligent data processing capabilities, Elegant 3 is easy to maneuver. Even beginners can get ideal scan results in minutes.




New UI

Neater and more user-friendly interface More flexible workflow

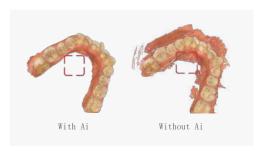

Al Scan

Automatically identifies and filters out unnecessary soft tissue data during scanning resulting in a quicker and cleaner scan process.

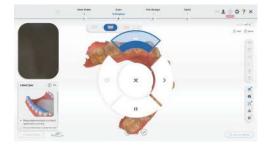

Motion Sensing


Dentists can scan with minimal contact with the computer, ensuring a safe and efficient scan experience.

Clinical Toolkit

The clinical toolkit assists dentists to evaluate and pre-design the scanned data in their own clinic. It facilitates effective communication with dental labs to deliver quality dental work with high efficiency.



Undercut Check

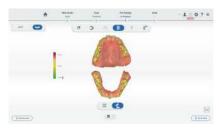

Bite Check



Margin Line Auto Extraction




Tooth Marking



Coordinate Adjustment



Direct Export into exocad

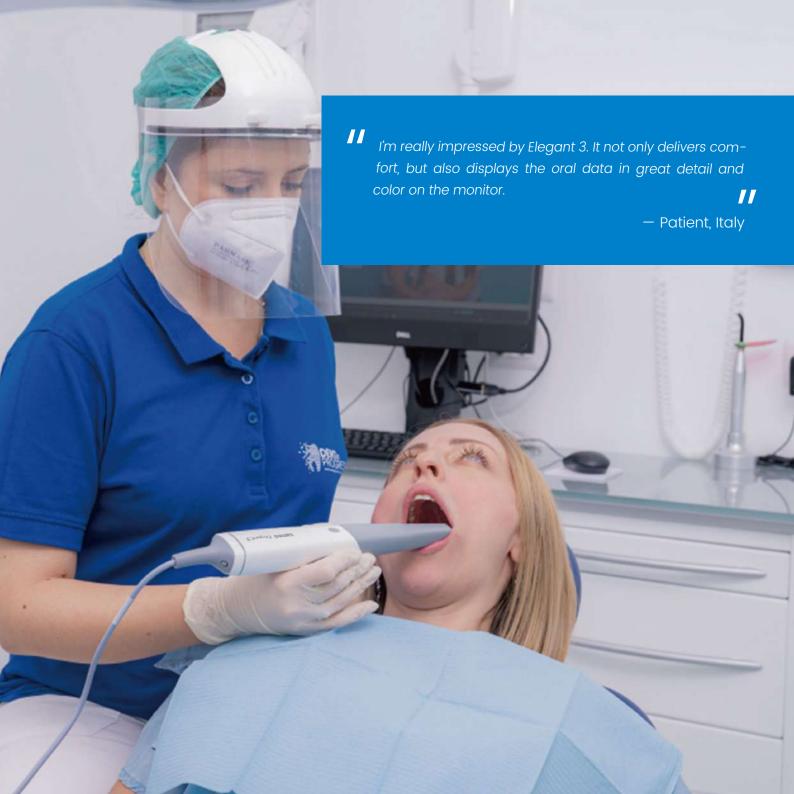

Intuitive Dentist-Patient Communication

Clear and Engaged – Satisfactory Experience for Patients

Elegant 3, capable of more than just digital impression acquisition, is integrated with powerful 3D visualization tools to help dentists engage patients with their treatment plan and improve overall satisfaction.



Oral Health Report

Health report includes information such as existing dental caries, missing teeth, dental calculus, pigment, etc. Patients have a better understanding of their oral health.




Ortho Simulation

Patients can preview images of their straightened-up results in real time.



Scan-to-print Process

Dental models can be edited and printed chairside enabling the direct visualization of the treatment plan and shortening the restoration fabrication cycle.

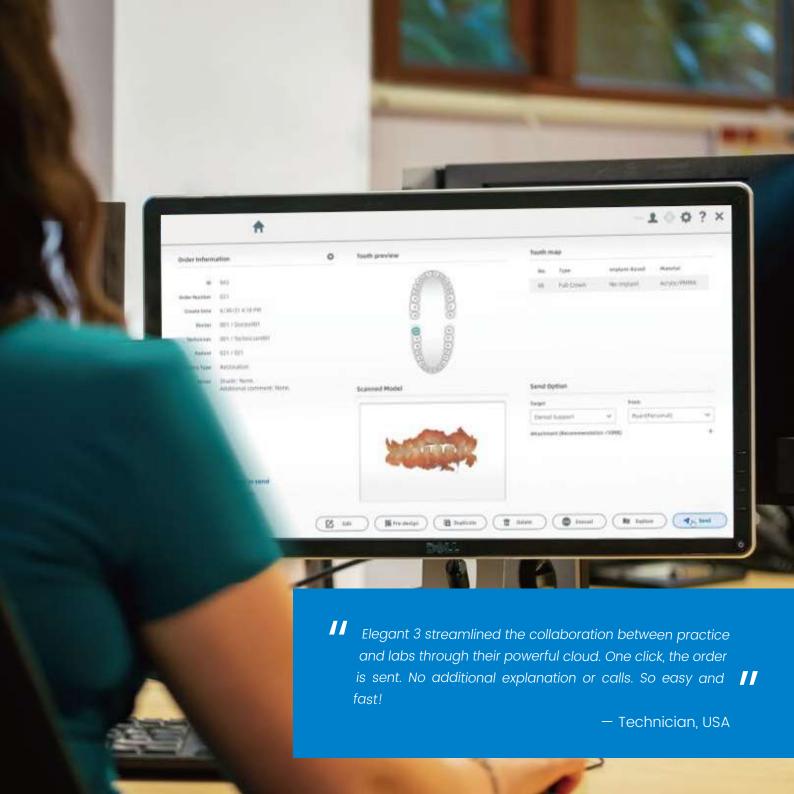

Efficient Dentist-Technician Collaboration

Standard and Unlimited - Seamless Collaboration for Dentists and Technicians

Elegant 3 helps streamline the collaboration between clinics and labs to achieve efficient and effective restorative works.

Cloud Data Platform

Cloud data platform enables data transmission, making it convenient for data transfer and communication between clinics and dental labs.



Extensive Range of Applications

Restoration, Implant, Orthodontics – Dental 3D Scanning Solutions

Elegant 3 intraoral scanner can be applied for clinical indications range spans general restorations, implants and orthodontics. Whether through same-day chairside restorations or clinic-dental lab workflow, it always provides a state-of-the-art user experience.



Orthodontics



Veneer



Inlay-onlay



6-Unit Bridge



All-on-6 Implant



Pediatric Dental Care





COOPERATIVE BUSINESS





















hyperDENT 3shape ▶

Technical Specifications

cameo Elegant 3

Scan Field	Standard scanner tip: 16 mm × 12 mm × 22 mm Mini scanner tip: 12 mm × 9 mm × 22 mm
Scan Depth	-2-20 mm from exit surface of tip
Scan Principle	Non-contact scanner with structured light
Dimension (L × W × H)	281 mm × 33 mm × 46 mm
Weight	240 ± 10 g (without cables)
Output	STL, OBJ, PLY
Connector	USB 3.0
Power Input	12 V DC/3 A

Recommended PC Configuration

Intel Core i7-8700 or higher
16 GB or higher
256 GB SSD or above
NVIDIA® RTX 2060 6GB or higher
Windows 10 Professional (64-bit) or later versions of Windows operating system
1920 × 1080, 60Hz or higher
More than 2 type-A USB 3.0 (or higher) ports



cameo Elegant 3



Aidite (Qinhuangdao) Technology Co.,Ltd.



Tel: 0086-335-8587898 Fax: 0086-335-8587198 Email: info@aidite.com

Email: info@aidite.com NMPA ISO13485 **C€**0197 **FDA**K111291

